

# After a LENA Day


## Process LENA devices with LENA Hub

To get help, visit [learn.lena.org/library/process](https://learn.lena.org/library/process).



### 1. Open LENA Hub.

If prompted, log on with your LENA username and password.

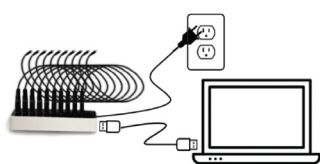

### 2. Set up one Multi Dock.

- Start with one empty Multi Dock and cords — no LENA devices attached.
- Connect the Multi Dock to a power outlet.
- Connect the Multi Dock to the LENA Hub.



### 3. Turn on all LENA devices and check screen message.

Device Message	What it Means	What to Do
<b>Paused</b> <b>xx:xx:xx</b> or <b>Memory Full</b>	Device contains data.	Process in Step 4.
<b>Paused</b> (no time shown)	Device is empty.	Set aside for next LENA Day.




### 4. Connect, confirm, repeat.

- Connect ONE device with data to the Multi Dock.
- Confirm the name or serial number in Hub.
- Repeat, one device at a time, until the Multi Dock is full.



### 5. Let connected devices transfer to Hub.

- Wait for all devices to show a battery icon.
- Disconnect the full Multi Dock from the computer, and leave it plugged into the power outlet to charge devices.
- Repeat steps 2-5 with another Multi Dock as needed.

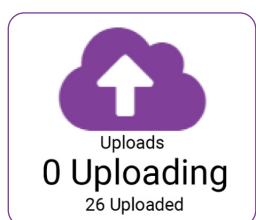

### 6. Leave Hub open until all LENA files upload to the cloud.

- Computer must be on, laptop lid open, connected to internet, and not sleeping until Hub shows 0 Uploading.
- See a gear icon instead of a cloud? Wait for Hub to show 0 Processing.